Wisbi®

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User manual

Mode d'emploi

Wisbi Home - WW1901

www.TEXIBLE.com

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SAFETY INFORMATION

Please carefully read through these instructions before using the device! TEXIBLE GmbH assumes no liability for material damages and personal injuries that arise due to non-observance of the operating instructions and especially the safety information!

TEXIBLE Wisbi® is not a medical device and must not be used to diagnose or treat illnesses. Ask your doctor for prevention, diagnosis and treatment of illnesses.

The TEXIBLE Wisbi transmitter must be placed under or next to the bed in a way that there is no danger for the patient through the cable.

Note for care workers

TEXIBLE Wisbi is a support aid for care workers and does not serve as a replacement. TEXIBLE Wisbi must be checked for functionality before each use. A constant internet connection is required.

Maintenance and troubleshooting

TEXIBLE Wisbi is free of maintenance except when the battery needs replacement. Only tighten the screws on the battery compartment slightly.

The following steps should be considered if the function is restricted:

- Is the bed pad correctly connected and the transmitter switched on?
- Is the TEXIBLE Wisbi transmitter connected with the WLAN?
- Is the transmitter battery functioning?
- Is the bed pad free of defects?

If the mentioned points are fulfilled, but the function is still limited, please contact the technical staff.

Note: The TEXIBLE Wisbi transmitter and the included QR code form an encrypted paging system. Keep the QR code secure and protect it from misuse. This ensures that only you have access to the information of your TEXIBLE Wisbi system.

Make sure that the radio link is unobstructed. Do not mount in distribution boxes, metal enclosures or close to large metal objects.

ATTENTION: The TEXIBLE Wisbi App cannot access the basic settings of your mobile device. Note, therefore, that no audible alarm sounds in the silent mode of the mobile device. Check the sound settings of the mobile device before each use.

The mobile device requires push notifications. You can test this by putting your mobile device into standby mode and performing a functional test.

USER MANUAL

1. PACKAGE CONTENTS

The package contents:

- 1x TEXIBLE Wisbi®-Bed Pad
- 1x TEXIBLE Wisbi®-Transmitter (2,4GHz 802.11b/g/n, WPA/WPA2 PSK)
- 1x Unique QR-Code for TEXIBLE Wisbi App

2. SCOPE OF APPLICATION AND FUNCTION

TEXIBLE Wisbi is the innovative nursing care product that stands out for its ability to automatically signal an alarm when moisture, body fluids or a bed exit are detected.

The sensor textile is used in a wide variety of nursing applications, such as senior and nursing homes, assisted living and family care. The necessary basic care, as well as the incontinence care and the supervision of demented and fall-endangered residents represents an increased care effort for all persons involved.

TEXIBLE Wisbi offers the right support here. The system allows an optimal skin climate and minimizes the risk of running away.

3. PUTTING TEXIBLE WISBI INTO OPERATION

The TEXIBLE Wisbi transmitter is equipped with a WLAN module WPA/WPA2 PSK (2,4GHz 802.11b/g/n) and is integrated into your home network.

Together with the TEXIBLE Wisbi App on your mobile device you have your personal call system. Every 4 seconds a measurement is carried out.

When detecting a wet bed pad or a bed exit, an alarm is sent to your mobile device via WLAN and an active Internet connection and you receive an acoustic alarm.

1. DOWNLOAD THE APP:

Download the App „TEXIBLE Wisbi“ from the Google Play or the App Store. Start the App and follow the manual.

2. **SCAN THE QR CODE:** Scan the QR code on the back of the TEXIBLE Wisbi transmitter with the TEXIBLE Wisbi App to pair the transmitter with your app


3. **INSERT BATTERIES:** Insert three mignon batteries (AA) into the battery compartment. When the batteries are inserted correctly, the white LED flashes briefly. Switch the TEXIBLE Wisbi transmitter on through pressing the button  once. The LED flashes whit briefly.

4. **CONNECTION BED PAD:** Place the bed pad as the top layer on the sheet. White cotton side up, blue side down. Then fasten the TEXTIBLE Wisbi transmitter using the push buttons with the bed pad. (See image below)

5. **WLAN PAIRING:** There are two ways to connect your Wisbi transmitter with your WLAN:

1 WITH THE WPS/WIFI BUTTON ON YOUR ROUTER:

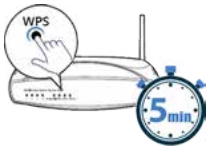Press the button  4 times within 5 seconds (LED flashes white) and then within 5 minutes press the WPS button on your WLAN router.

The TEXTIBLE Wisbi Sender now automatically connects to your home network. If the connection is OK, the white LED glow for 3 seconds.

After successful pairing, it can take up to 30 seconds until the status message in the app is updated.

In the future, the transmitter will automatically connect to your home network after you turn it off and on again.

2 ENTER YOUR WLAN NAME AND PASSWORD MANUALLY

Press the key  6 times within 5 seconds (LED flashes red).

Open the WLAN settings on your mobile device and connect to WLAN „Wisbi“.

Now enter within 5 minutes the information (SSID / WLAN name and Password) of your home network in the app.

The TEXTIBLE Wisbi Sender now automatically connects to your home network. If the manual connection is OK the red LED glow continuously for 3 seconds.

6. **SENSOR SETTINGS:** Configure the settings of the system, such as the delayed alert when leaving the bed, to suit your individual care needs.

4. FUNCTIONAL TEST

TEXTIBLE Wisbi must be checked for proper functioning before each use. Switch on the TEXTIBLE Wisbi transmitter. Press the button  for 3 seconds until the orange LED flashes. This will lead to a functional test of the electronics and sends an alarm.

5. ACKNOWLEDGE ALARMS

When an alarm is triggered, a status message appears in the drop-down menu on your mobile device. To confirm this alarm you have to tap on this status message. Otherwise the alarm repeats every 30 seconds.

6. SWITCH OFF

Once TEXTIBLE Wisbi is switched on it is constantly in standby mode and measures every 4 seconds. If the transmitter is not used, it can be switched off by pressing the button  for 3 seconds to preserve the batteries. The LED glows red. By pressing the button  again once, the transmitter is switched on again.

7. OPERATION AND SETTINGS APP


If your supervised person is lying on the TEXTIBLE Wisbi bed pad for at least 10 seconds, occupancy detection will be activated.

In the settings of the App you can also set a delay up to 15min. For example, your supervised person can independently go to the bathroom. An alarm will be sent only after this time runs up and the person did not go back to bed.

You can enable and silence alarms. Important is, that your mobile device is not in silent mode.

8. CHANGING BATTERIES OF TRANSMITTER

Open the rear cover of the battery compartment by loosening the screw. Insert three mignon batteries (AA) into the battery compartment as shown and close the cover again.

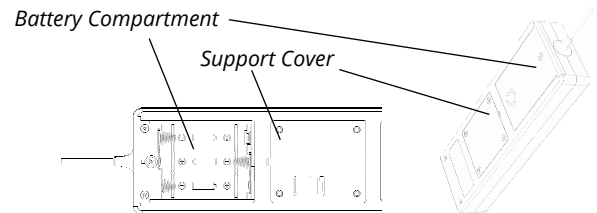9. OPERATING CONCEPT

Function	Button	LED Display
Switching TEXIBLE Wisbi Transmitter ON	 Press 1 x	flashing white
Switching TEXIBLE Wisbi Transmitter OFF	 Press 3 sec	glowing red
Functional test	 Press 3 sec	flashing orange
Linking WPS WLAN	 Press 4 x	flashing white (connected: constant light)
Manual WLAN setup	 Press 6 x	flashing red (connected: constant light)
Alarm		flashing orange
Battery low		red & orange flash alternately
Error		glowing red

10. DISPOSAL INFORMATION

Dispose the TEXIBLE Wisbi transmitter with accompanying batteries at a designated collection point.

The TEXIBLE Wisbi bed pad can be disposed in household waste. Please contact your waste disposal company or the representative for more information about nearby collection points.

Dispose of the packaging material in the collection containers for cardboard, paper and plastics.

11. WARRANTY

We rectify any device defects due to material or manufacturing faults by repair or replacement free of charge within the statutory warranty period. The warranty expires in the case of intervention by third parties.

12. EC - DECLARATION OF CONFORMITY

Manufacturer:
TEXIBLE GmbH, Schwefelbadstrasse 2, 6845 Hohenems, Austria

The following guidelines have been complied with:
EMC 2014/30/EU, RED 2014/53 EU, LVD 2014/35/EU,
RoHS 2011/65/EU

The declaration of conformity can be requested from the manufacturer.

13. CUSTOMER SERVICE

Please contact the manufacturer or your dealer if malfunctions still occur despite proper handling or if the device has been damaged.

14. TECHNICAL DATA TEXIBLE WISBI

Size	64 x 100cm
Weight	250g
Water absorption	up to 700ml/m ²
Trigger quantity	10ml
Minimal body weight (for occupancy detection)	45kg
Material	67% Cotton 30% Polyester 3% Polyurethane
Wash resistance	    
Power supply	3xAA mignon batteries, LR6 
Power	1,3mW
Protection class	IP40
Ambient temperature	0 to +55°C
Battery life	1 year depending on application
WLAN Connection	2,4GHz 802.11b/g/n WPA/WPA2 PSK

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TEXIBLE GmbH, Schwefelbadstrasse 2, 6845 Hohenems, Autriche